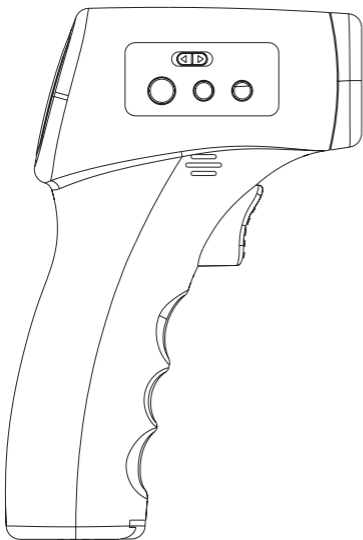
© copyright Edco Deutschland GmbH
Siemensstr. 31 47533 Kleve
Deutschland. Made in P.R.C.
www.grundigaccessories.com
© The Grundig trademark is used
under license

Infrared Thermometer

Instruction Manual

Customer Model:MDI907 Factory Model:MDI903

Tip: A wrong operation can cause a inaccurate measurement result, the correct measured distance is 1-5 cm to forehead (between the eyebrows).

July 2020 VER.1.0

Table of Contents

Introduction	2
Safety Information	3
Product structure	6
Battery installation	7
Setting	8
Use of the unit	9
Trouble shoot	12
Care and maintenance	13
Specification	14
Warranty information	15

Introduction

▲ Your new non-contact forehead thermometer uses advanced infrared (IR) technology to measure temperature instantly and accurately on the forehead or object.

▲ Easy to use and less measurement time. This thermometer does not need to contact body or object to ensure the safety and hygiene. Its ergonomic design makes this thermometer be simple and very easy to check the temperature. It only takes 1 second to take measurement and reading.

▲ Body and object modes

This thermometer supports to measure temperature of body and object. Measuring range of object mode is 0.0 to 100.0 °C(32.0 °F~212.0 °F). That means except to take body temperature, it also supports to take the surface temperature of object as below,

.Surface temperature of milk in the bottle

.Surface temperature of baby bath

.Environment temperature

▲ Alarm sounds when body temperature is over 38.5 °C, an alarm sounds for 3 times.

▲ Memory function

It stores up to 1 sets recent measurement data.

▲ Auto power off

Automatically power off if left idle for 30 seconds .

▲ Measurement

Measurement time interval in 15s and the measurement distance is within 1-5 cm.

Please read the manual carefully before you use the unit, and keep for future reference.

⚠ This thermometer is not intended to substitute for a consultation with your physician. The forehead scan temperature serves as a reference only. It cannot be a judgment on fever.

⚠ Basic safety precautions should always be observed, especially when the thermometer is used on or near children and disabled persons.

⚠ Please place the device unreachable by young.

⚠ Avoid direct sunlight.

⚠ Do not touch the lens.

⚠ No modification of this device is allowed.

⚠ The swallowing of small parts like packing bag, battery, battery cover and so on may cause the suffocation.

⚠ Please do not use a dilution agent, alcohol or petrol to clean the unit. Please treat it gently and prevent the falling from a high place.

⚠ Please do not immerse it in liquid.

⚠ Never leave battery in the battery compartment for a long time without use, as they may leak and cause damage to the unit.

 Please take off the battery if you do not intend to use within 3 months. Replace with new batteries if the unit display a low battery symbol.

 Do not mix the old and new batteries together.

 This article is not for professional usage.

Safety Information

▲ To ensure the correct use of the product, basic safety measures should always be followed including the warning and the caution symbols listed in the instruction manual :

Symbol descriptions	
The following symbols may appear in this manual, on the label, on the device, or on its accessories. Some of the symbols represent standards and compliance associated with the device and its use.	
	WARNING :This alert identifies hazards that may cause serious personal injury or death .
	CAUTION :This alert identifies hazards that may cause minor personal injury, product damage, or property damage .
	Type BF applied part
	Manufacturer
SN	Specifies serial numbe
	DISPOSAL : Do not dispose this product as unsorted municipal waste. Collection of such waste separately for special treatment is necessary.
	Direct current
	Operating instructions
	Consult accompanying documents

Safety Information

WARNING :

Do not dispose of electrical appliances as unsorted municipal waste, use separate collection facilities. Contact you local government for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being.

Classification

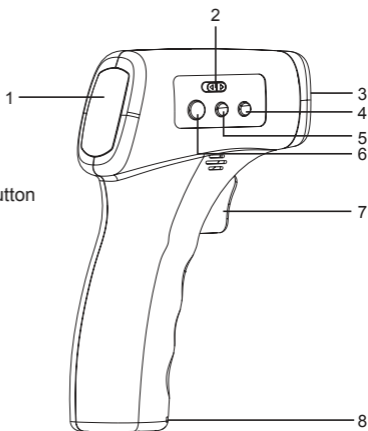
- 1.Internally powered equipment;
- 2.Type BF applied part;
- 3.Protection against ingress of water: IPX0;
- 4.Not category AP/APG equipment ;
- 5.Mode of operation : Continuous operation ;

Note : the user must check that the equipment functions safely and ensure that it is in proper working condition before it being used.

Product structure


Body

1. Display screen
2. Mode switch
3. Sensor probe
4. Down (-) button
5. Up (+) button
6. Setting button
7. START/STOP button
8. Battery cover

Display

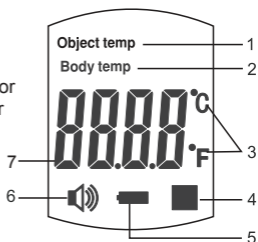1. Object temperature indicator
2. Body temperature indicator
3. Temperature unit
4. Memory mode
5. Low battery indicator
6. Sound indicator
7. Temperature reading

Battery Installation

1. Remove the battery cover from the battery compartment as the arrow direction accordingly.
2. Insert 2 pcs 1.5 V AAA size alkaline batteries, and ensure each battery is in the proper direction. Close the battery cover.

Low battery and replacement

During power on, the low battery symbol  will display once the unit starts to work, and you should replace with new batteries, otherwise the unit can not work.

Battery type and replacement

Please use 2 pcs of AAA identical 1.5V alkaline batteries. Do not use the batteries beyond their expiry date .

Please remove the batteries if you do not need to use for long time .

WARNING

Dispose the battery in accordance with all federal, state and local laws .To avoid fire and explosion hazard, do not burn or incinerate the battery.

How to set

Display	Function	Description
F-1	Unit selection	Toggle between Unit °C and °F
F-2	Turn On/Off the alarm sound	With "OPEN" on the screen means the sound is on With "CLOS" on the screen means the sound is off

Use of the unit

Pre-measurement

About normal body temperature & fever the temperature in the forehead and temple area differs from the internal temperature, which is taken orally or rectally. Vasoconstriction, an effect which constricts the blood vessels and cools the skin, can occur during the early stages of a fever. In this case, the temperature measured by the non-contact thermometer may be unusually low. If the measurement therefore does not match the patients own perception or is unusually low, repeat the measurement every 15 minutes. As a reference, you can also measure the internal body temperature using a conventional oral or rectal thermometer.

Body temperature can vary from one individual/person to next. It also varies by location on the body and time of day. Below shows the statistical normal ranges from different sites.

Please keep in mind that temperatures measured from different sites, even at the same time, should not be directly compared. Fever indicates that the body temperature is higher than normal. This symptom may be caused by infection, overdressing or immunization. Some people may not experience fever even when they are ill.

These include, but are not limited to, infants younger than 3 months old, persons with compromised immune systems, persons taking antibiotics, steroids or antipyretics (aspirin, ibuprofen, acetaminophen), or persons with certain chronic illnesses. Please consult your physician when you feel ill even if you do not have fever.

Table* 1 Body Site Normal Temperature Range


Body site	Normal temperature range
Oral	0.6°C (1°F) or more
Rectal/ear	0.3°C to 0.6 °C (0.5 °F to 1°F) higher than oral temperature
Axillary (armpit)	0.3°C to 0.6°C(0.5°F to 1°F) lower than oral temperature

Use of the unit

Measuring procedure

As a body thermometer

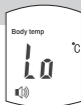1) Select body mode , the temperature unit will appear. (Picture at right)

2) Move the probe close to the forehead, make sure the probe is flat and take measurement with a distance within 1-5 cm. Press Start/Stop button.

3) The measurement will be finished in 1 seconds, when it has been completed, you can hear 3 short beeps, at the same time, the reading will appear on the display.

Note :

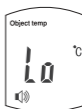
- If the reading is $<37.5^{\circ}\text{C}$ (99.5°F) and $\geq 32^{\circ}\text{C}$ (89.6°F) the display will show with green backlight.
- If the reading is between value $37.5\text{--}38.5^{\circ}\text{C}$ ($99.5\text{--}101.3^{\circ}\text{F}$), the display will show together with 1 short beeps, the orange backlight lasts for alarm.
- If the reading is \geq above value 38.5°C (109.4°F), the display will show together with 3 short beeps, the red backlight lasts for alarm.
- As the forehead measurement temperature is likely to be affected by sweat, oil and the surroundings, the reading shall be taken as a reference only.
- If the probe is placed at an angle close to the forehead measurement, the reading will be affected by surrounding temperature.
- Babies' skin reacts very quickly in the ambient temperature. Therefore, do not take their temperature with the non-contact thermometer during/after breastfeeding, because the skin temperature maybe lower than the internal body temperature.
- The thermometer will automatically turn off if left idle for 30 seconds.

Use of the unit

Measuring procedure

As a object thermometer

1. Select Object Mode, the temperature unit will appear. (Picture at right)
2. Move the probe close to the object, make sure the probe is flat and take measurement with a distance within 1-5cm. Press the Start/Stop button to take measurement
3. The measurement will be finished in 1 seconds, when it has been completed, you can hear 3 short beeps, at the same time, the reading will appear on the display






Note :

- As the object scan temperature is likely to be affected by the surroundings, the reading shall be taken as a reference only.
- If the probe is placed at an angle close the object scan, the reading will be affected by surrounding temperature.
- The thermometer will automatically turn off if left idle for 30 seconds.

Trouble shoot

Error indicators

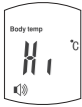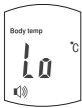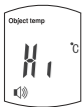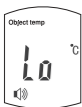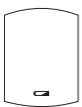The following symbols will appear on the display when measurement is abnormal

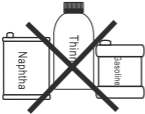
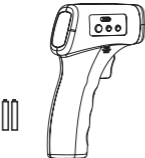
Symbol	Correction
	In body mode. Measured temperature is above measuring range of $>43.0^{\circ}\text{C}$ (109.4°F) It is shown on the LCD with red LED.
	In body mode. Measured temperature is below Measuring Range of $<32^{\circ}\text{C}$ / $89,6^{\circ}\text{F}$ It is shown on the LCD with red LED.
	In object mode. Measured temperature is above measuring range of $>100^{\circ}\text{C}$ / 212°F or environmental temperature is above system operating range $>40^{\circ}\text{C}$ / 104°F It is shown on the LCD with red LED.
	In object mode. Measured temperature is below measuring range of $<0.0^{\circ}\text{C}$ / $32,0^{\circ}\text{F}$ Or environmental temperature is below system operating range $<5^{\circ}\text{C}$ / 41°F It is shown on the LCD with red LED.
	If low battery indicator appears, replace all batteries with new ones.
ERR	Thermomter sytem fails or affected by electric magnetic field.

Care and maintenance

Care for the main unit

- Keep the unit in the storage case when not in use.
- Clean the unit with a soft dry cloth. Do not use any abrasive or volatile cleaners.
- Never immerse the unit or any of its component in water.

Maintenance

<p>Do not clean the unit with naphtha, thinner or gasoline etc.</p>	<p>Store the unit in a clean and dry location. Do not subject the unit to extreme hot or cold temperature, humidity and direct sunlight.</p>
	
<p>Remove the batteries if the unit will not be used in 3 months or longer.</p>	<p>Do not use the unit under the influence of electromagnetic interference (nearby cellphones, microwave etc.)</p>
	

Note : We will not be responsible for any quality problem if you do not care and maintain the product as instructed.

Specification

Description	Non-contact Infrared Thermometer
Display	LCD digital display
Measuring localization	Forehead and object surface
Measurement range	Body mode 32.0°C~43.0°C(89.6°F~109.4°F); Object mode 0.0~100.0°C(32.0°F-212.0°F);
Temperature unit	°C/°F
Display-resolution	0.1°C/0.1°F
Accuracy	±0.2°C/±0.4°F (within 35.0°C~42.0°C/95.0°F~107.6°F)
Memory function	Last memory of measurement values
Power source	2 pcs AAA alkaline battery
Automatic power-off	In 30 seconds
Main unit weight	Approx. 93 gr (batteries not included)
Main unit size	95mm*45mm*155mm
Battery life	Could be used for 500 times for normal condition
Accessories	Instruction manual + storage pouch
Operating environment	Body mode 10.0°C~40.0°C/50.0 °F -104.0°F Object mode 5°C~40.0°C/41.0°F~104.0°F
Storage environment	Temperature >-20°C~+55°C/-4.0°F ~+131.0°F Humidity: RH≤ 85% does not condense avoid crash, sunburn or rain during transportation

Warranty information

▲ The unit is guaranteed to be free of defects in workmanship and materials under normal use for a period of 1 year from the date of purchase.

▲ For repair under this warranty. Our authorized service agent must be advised of the fault within the period of the warranty. This warranty only covers parts and labor service under normal operations. Any defect resulting from natural causes, eg. flood, hurricane etc, is not covered in this guarantee. This guarantee also does not cover damage incurred by use of the unit not in accordance with the instructions, accidental damage, or being tampered or serviced by unauthorized service agents.

▲ The following will be excluded from this warranty-if the thermometer has been misused, abused, or neglect in following the manual's instructions on purpose and unauthorized repair or modifications.

 The device requires no calibration.

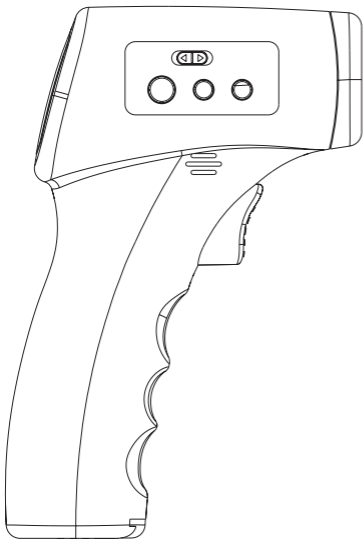
 The device is not repairable and contains no user serviceable parts

 ShenZhen Medek Bio-Medical Co.,Ltd.
ADD: No. 101 Wenhao Industrial Park No.13 Yuanxin Road
Tongle Community, Baolong Street, Longgang District. 518116
Shenzhen P.R.China

 Name:Shanghai International Holding Corp.GmbH(Europe)
Add: Eiffestrasse 80,20537 Hamburg,Germany

© copyright Edco Deutschland GmbH
Siemensstr. 31 47533 Kleve
Deutschland. Made in P.R.C.
www.grundigaccessories.com
© The Grundig trademark is used
under license

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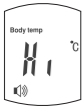
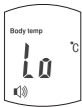
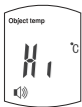
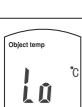
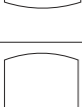
Symbol	Działanie
 The image shows a rectangular LCD display with a white background. At the top left, it says "Body temp". In the center, the letters "Hi" are displayed in a large, bold font. To the right of "Hi" is a small "°C" symbol. At the bottom left, there is a small speaker icon with sound waves. A red LED indicator is shown as a small red rectangle at the bottom right of the display.	<p>W trybie ciało. Mierzona temperatura przekracza zakres pomiaru $>43.0^{\circ}\text{C}$(109.4°F) Wyświetlacz LCD zostaje podświetlony przez czerwoną diodę LED.</p>
 The image shows a rectangular LCD display with a white background. At the top left, it says "Body temp". In the center, the letters "Lo" are displayed in a large, bold font. To the right of "Lo" is a small "°C" symbol. At the bottom left, there is a small speaker icon with sound waves. A red LED indicator is shown as a small red rectangle at the bottom right of the display.	<p>W trybie ciało. Mierzona temperatura jest poniżej zakresu pomiaru $< 32^{\circ}\text{C}$/$89,6^{\circ}\text{F}$ Wyświetlacz LCD zostaje podświetlony przez czerwoną diodę LED.</p>
 The image shows a rectangular LCD display with a white background. At the top left, it says "Object temp". In the center, the letters "Hi" are displayed in a large, bold font. To the right of "Hi" is a small "°C" symbol. At the bottom left, there is a small speaker icon with sound waves. A red LED indicator is shown as a small red rectangle at the bottom right of the display.	<p>W trybie przedmiot. Mierzona temperatura przekracza zakres pomiaru $> 100^{\circ}\text{C}$/212°F Lub temperatura otoczenia jest powyżej systemowego zakresu działania $> 40^{\circ}\text{C}$/104°F Wyświetlacz LCD zostaje podświetlony przez czerwoną diodę LED.</p>
 The image shows a rectangular LCD display with a white background. At the top left, it says "Object temp". In the center, the letters "Lo" are displayed in a large, bold font. To the right of "Lo" is a small "°C" symbol. At the bottom left, there is a small speaker icon with sound waves. A red LED indicator is shown as a small red rectangle at the bottom right of the display.	<p>W trybie przedmiot. Mierzona temperatura jest poniżej zakresu pomiaru $< 0,0^{\circ}\text{C}$/$32,0^{\circ}\text{F}$ Lub temperatura otoczenia jest poniżej zakresu roboczego systemu $<5^{\circ}\text{C}$/41°F Wyświetlacz LCD zostaje podświetlony przez czerwoną diodę LED.</p>
 The image shows a rectangular LCD display with a white background. At the bottom center, there is a small icon of a battery with a low level of charge, represented by a short horizontal bar.	<p>Niski poziom baterii wymień zużyte baterie na nowe.</p>
ERR	<p>Błąd systemu lub zakłócenie działania przez pole elektromagnetyczne.</p>

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