Digitale BodyFit-Analysewaage
Digital BodyFit scale
Échelle d'analyse BodyFit numérique
Digitale BodyFit-analyseschaal
Cyfrowa skala analizy BodyFit
Escala digital de análisis BodyFit

BS-310.1FW

DE

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Digitale BodyFit-Analysewaage

BS-310.1FW

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
BS-310.1FW

Parts Description

1. Set button
2. Button
3. Button
4. Battery compartment (located at back of scale)
5. UNIT selection button (located at back of device)
6. Screen
7. Foot placement area

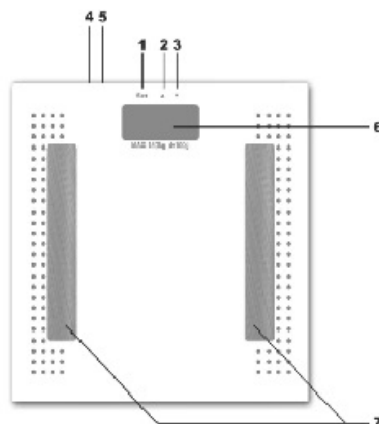Safety instructions:

Please read carefully before use

- Place the scale on a firm, flat and dry surface.
- Always position your dry, bare feet flat on the metallic part of the scale.
- Step onto the scale with dry, bare feet.
- Do not drop the scale as this could damage it.
- Clean the scale with a damp, soft cloth and wipe dry straight away.
- This scale is not suitable for business or professional use.
- The scale is not a playing. Do not allow children to play with the scale.
- This scale is not suitable for use by anyone fitted with an electrical device such as a pace maker or similar.
- Never use a combination of new and used batteries.
- Always respect the polarity.
- Do not clean the surface of the scale with abrasive substances or abrasive cloths.

Batteries

Batteries: 2 X AAA included.

Insert the batteries in the battery compartment (4) located at the back of the scale, observing the polarity. The display mode is active.

If "LO" appears on the screen, replace both batteries with new ones to restore the normal display.

Instructions for use

To change the unit of measure:

In weighing mode (weight display), press "UNIT" (5) located at the back of the scale to select the required unit of measurement.

To check your weight:

To initialise, press the top of the scale until the screen (6) displays "0". The scale is now ready for use.

To check your weight, stand on the scale, the display will flash. Stay on the scale until the display stabilizes to show your weight.
The scale will switch off automatically after a few seconds.

To store user data

- Press “Set” (1).
- Press “ ▼ ” or “ ▲ ” to select the desired memory.
- Confirm your choice by pressing “Set” (1). The symbol “male ” or “female” will flash.
- Press “ ▼ ” or “ ▲ ” to select gender.
- Confirm your choice by pressing “Set” (1) and “age” will flash.
- Press “ ▼ ” or “ ▲ ” to select age.
- Confirm your choice by pressing “Set” (1) and “height” will flash.
- Press “ ▼ ” or “ ▲ ” to select height.
- Confirm your choice by pressing “Set” (1) to display “0.0”

To check your weight and calculate the ratio of fats, fluids, muscles

- Press “Set” (1), the scale will light up.
- Select your memory by pressing the “ ▼ ” or “ ▲ ” arrows to display your personal data.
- Stand on the scale with dry, bare feet. Place each foot on the foot placement area (7) and stand still. Your weight will be displayed. Small cubes will scroll down the screen. As soon as the rates appear (percentage of body fat-“FAT”, percentage of body fluids-“WATER”, percentage of muscle-“MUSCLE”, percentage of bone mass-“BONE” and the number of calories your body needs to stay within these percentages), step off the scale and make a note of the different rates displayed (they will be shown three times).
The scale will switch off automatically.

Waste Management

The symbol above and on the product means that the product is classed as Electrical or Electronic equipment and should not be disposed with other household or commercial waste at the end of its useful life.

The Waste of Electrical and Electronic Equipment (WEEE) Directive (2012/19/EU) has been put in place to recycle products using best available recovery and recycling techniques to minimise the impact on the environment, treat any hazardous substances and avoid the increasing landfill. Contact local authorities for information on the correct disposal of Electrical or Electronic equipment.

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